

# QP-CryoSwitch Controller User Guide

## Preparation

### GIT repository

For user convenience, a small repository containing all the necessary files can be found [here](#).

Inside the repository, you will find a copy of the installation guide (this file), the controller datasheet, and finally a library containing python code to interface with the CryoSwitch Controller.

### Python setup

The controller's library doesn't need any specific python distribution, but some extra python libraries must be installed.

The easiest way to do so is by using 'pip'. Inside the repository folder and using your base or virtual environment run:

```
> pip install -r requirements.txt
```


### Preparation Summary

- Clone the GIT [repository](#).
- Install the necessary python dependencies.

## Connections and ports information


The QP-CryoSwitch Controller is specifically designed for driving inductive loads at room or cryogenic temperatures.

The controller supports a maximum of 4 ports (version dependent) with 7 pins per port. The following figure shows the back plate and pinout of the controller.

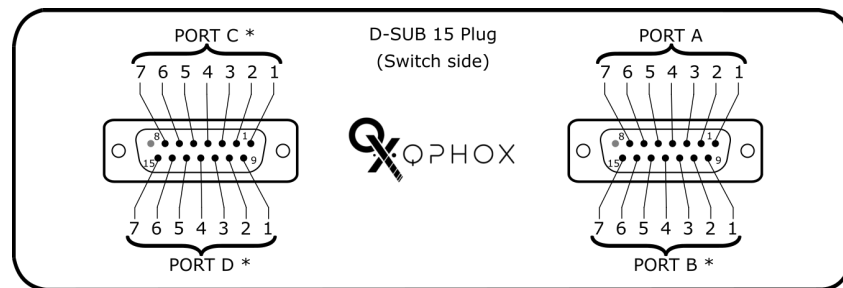

**Figure 1:** Controller output side

For example, if an SP6T cryogenic relay/actuator is used the connections will depend on the specific model. Two different connection examples are shown in the next figure:

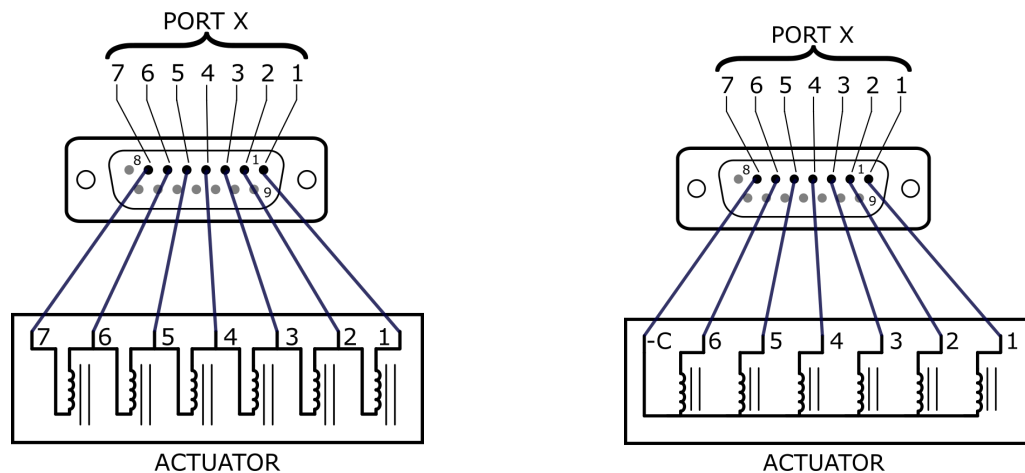

**Figure 2:** Application example

As seen in the first image of figure, 7 pins are used for driving a 6 contact actuator. Each pair of pins (1-2, 2-3..., 6-7) is used to drive a single relay contact of the actuator. In the second image, 6 contact pins and a common pin are used for driving the individual relay contacts.

## Python Library

The library contains the 'CryoSwitch' class which is a user-friendly way of interfacing between the QP-CryoSwitch Controller and a computer or another USB/Ethernet capable device.

The Cryoswitch class initialization takes some optional arguments:

**Cryoswitch(debug=False, COM\_port="", IP=None, SN=None)**

A basic implementation of the CryoSwitchController class can be done with the following functions:

**- start()**

- Input: None
- Default: None
- Enables the voltage rails, voltage converter and output stage

**- set\_output\_voltage(Vout)**

- Input: Desired output voltage (Vout)
- Default: 5V
- Sets the converter voltage to *Vout*. The output stage later utilizes the converter voltage to generate the positive/negative pulses.

**- set\_pulse\_duration\_ms(ms\_duration)**

- Input: Pulse width duration in milliseconds (ms\_duration).
- Default: 15ms.
- Sets the output pulse (positive/negative) duration in milliseconds.

**- connect(port, contact)**

- Input: Corresponding port and contact to be connected. *Port*={A, B, C, D}, *contact*={1,...,6}
- Default: None.
- Connects the specified switch contact of the specified controller port (switch).
- Returns: Current waveform.

**- disconnect(port, contact)**


- Input: Corresponding port and contact to be disconnected. *Port*={A, B, C, D}, *contact*={1,...,6}
- Default: None.
- Disconnects the specified switch contact of the specified controller port (switch).
- Returns: Current waveform.

## Typical application case

The CryoSwitch controller pulsing features use precise voltage and timing control to reduce heat dissipation, this way every user can optimize power dissipation while engaging the actuator.

Traditional electro-mechanical relays utilize an applied electromagnetic field in order to actuate a switch and offer high isolation and low cross-talk. These switches consist of an electromagnet in the form of a coil and a moving ferromagnetic plate that responds to the field in order to open or close the switch contact. The magnetic field is proportional to the coil drive current which must exceed a threshold in order to fully engage the actuator. Temperature differences (mK or room temp.), connectors, and cabling resistance will influence the equivalent line resistance ( $R_{line}$ ), hence limiting the current for a certain voltage pulse. **Figure 3** shows the equivalent diagram of a typical application including the line resistance. It's important to note that the line resistance  $R_{line}$  is not an actual resistor, but an equivalent for the entire system, thus a pulse configuration that works well for a room temperature test setup may not be optimal for a cryogenic setup where minimized power dissipation is critical.

By changing the output voltage and pulse width, it's possible to optimize the pulse current and compensate for the effective line resistance for either environment.

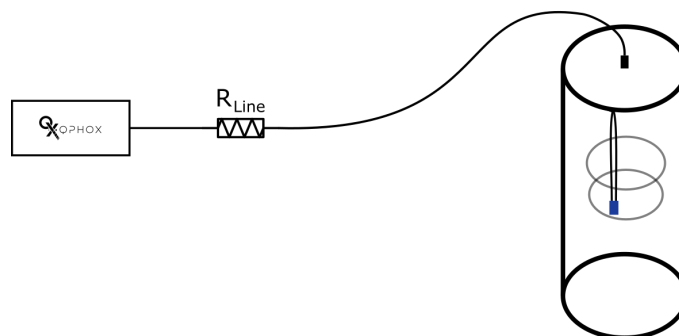

**Figure 3:** Typical cryogenic setup diagram

Going back to the [SDK](#), a typical application would look like this:

First, we initialize the controller software by declaring an instance of the Cryoswitch class:

```
switch = Cryoswitch()
```

There are two ways the user can connect to it, the first one is via ethernet. If ethernet is used, the user needs to provide the appropriate IP address. By default the IP address is 192.168.1.101 and the subnet mask 255.255.255.0, they can be easily changed by running the functions described in

if USB is used to establish the connection with the controller a COM port or serial number can be specified. Otherwise, the program will scan for available Cryoswitch controllers and connect to one.

Then we start the controller hardware by enabling the different voltage rails and setting the protections. By default, the output voltage is set to 5V and the pulse width to 15ms. This can be done by running the **start()** command:

**switch.start()**

At this point the controller is running, the last step before switching the actuator is setting the appropriate output voltage. As mentioned earlier, the appropriate pulse voltage depends on the user's setup configuration; however, as a starting point, we recommend a 5V pulse for cryogenic operation and a 20-25V pulse for room temperature operation.

An easy way to check if the actuator was engaged is by looking at the current waveform. To do this, the plotting function can be enabled by setting


**switch.plot = True**

The plotting function is enabled by default.

**Figure 4** shows a typical current waveform at room temperature.

When the pulse starts, the coil is energized and the rising current exerts a force on the actuator. When the force on the actuator exceeds the retaining force, the actuator starts to move. The motion of the actuator (usually made from a ferromagnetic material) induces a back electromagnetic force in the coil, effectively reducing the current. As the movement continues, the current continues to drop until the actuator reaches the final position.

Once the actuator is at rest, the current starts to rise, saturating the coil until it reaches its final DC current (given by the equivalent line resistance and pulse voltage). This part of the switching is not strictly necessary, therefore it's called the excess pulse. Although one might try to reduce the pulse width to avoid the excess pulse and therefore further power dissipation, it's difficult to know beforehand what the pulse duration should be. As discussed earlier, this depends on the fridge wiring and temperature. Having said that, pulse widths of 15-20ms are usually enough to engage the actuator.

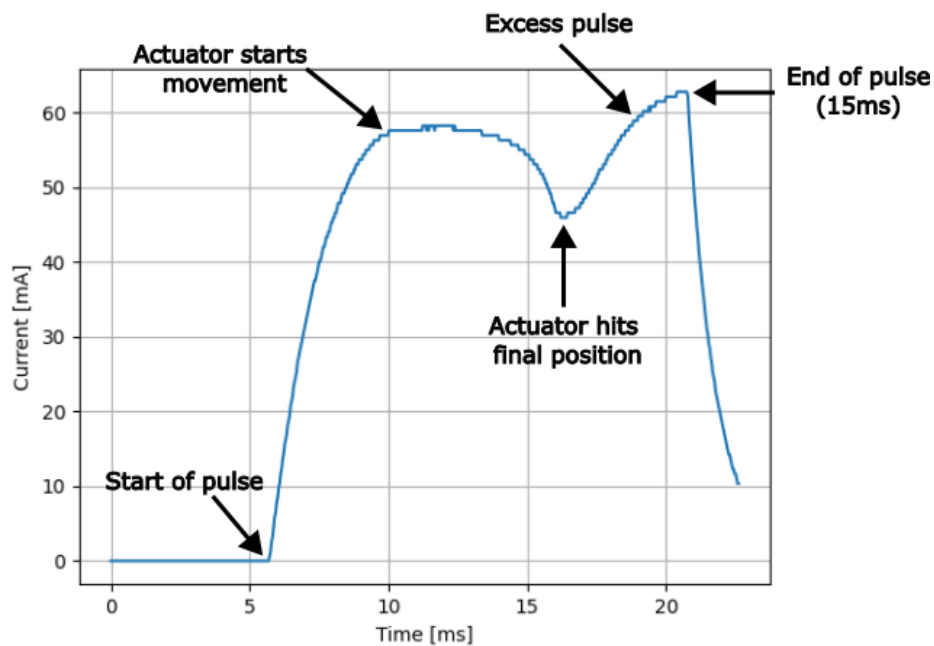


**Figure 4:** Normal operation waveform

## Troubleshooting:

### Recognizing insufficient output voltage:


The traces in **Figure 5** show a fully engaged actuator waveform (**blue**) and a non-engaged actuator waveform (**orange**). Since the current is directly proportional to the pulse voltage, reducing the pulse voltage will also reduce the pulse current and vice versa. The actuator needs a minimum amount of current to overcome the actuator retaining force (around 60mA). A pulse that does not provide sufficient current to overcome the actuator retaining force will show a waveform characteristic of the **orange** trace. This can be simply solved by increasing the pulse voltage.

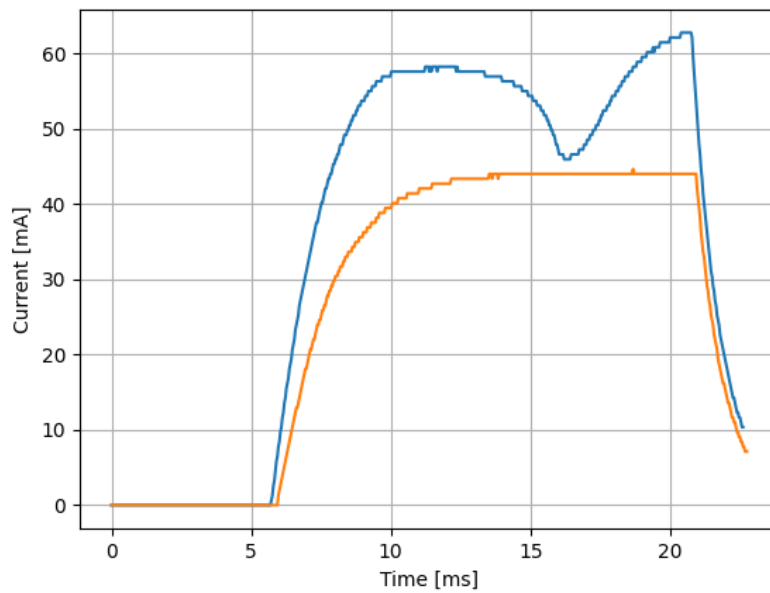

**Figure 5:** Low voltage waveform

### Recognizing open-circuits:

Another common scenario is unconnected/broken cables, either inside the fridge or along the line. Internally, for every pulse, the controller uses some of the pulse current to bias the internal output stage. The bias current depends on the pulse voltage and is between 5-25mA. **Figure 6** shows a typical waveform when the outputs are not connected. Although the plot shows there is some current flowing (bias current), the user can easily distinguish this from an actual pulse and infer that the output is probably not properly connected.

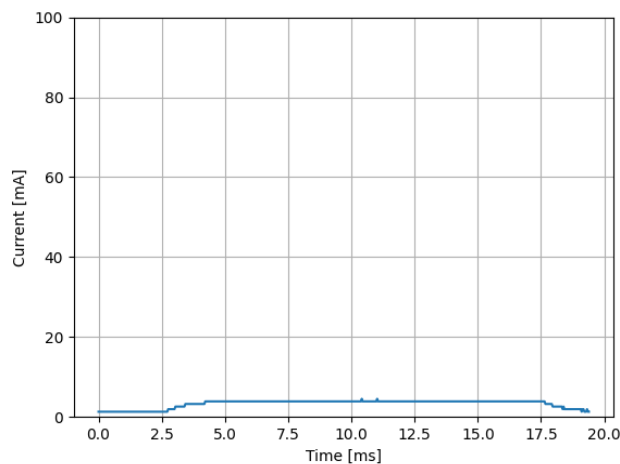

**Figure 6:** Unconnected waveform

## Recommendations:

In summary, the user should set the appropriate pulse voltage and pulse width depending on their fridge configuration but a good starting point would be:

Setup	Voltage	Width
Room temperature	20V	15ms
Cryogenic	5V	10ms
Recommended settings		

The command sequence discussed above can be summarized as:

```
switch = Cryoswitch() ## -> CryoSwitch class declaration and USB connection
switch.start() ## -> Initialization of the internal hardware
switch.plot = True ## -> Enable the current plotting function
switch.set_output_voltage(5) ## -> Set the output pulse voltage to 5V
```



## Overcurrent related functions

### - enable\_OCP()

- Input: None
- Default: None.
- Enables the overcurrent protection.

### - set\_OCP\_mA(OCP\_value)

- Input: Overcurrent protection trigger value (OCP\_value).
- Default: 100mA.
- Sets the overcurrent protection to the specified value.

### - enable\_chopping()

- Input: None.
- Default: None.
- Enables the chopping function. When an overcurrent condition occurs, the controller will 'chop' the excess current instead of disabling the output. Please refer to the installation guide for further information.

### - disable\_chopping()

- Input: None.
- Default: None.
- Disables the chopping function. When an overcurrent condition occurs, the controller will disable the output voltage. Please refer to the installation guide for further information.

### - reset\_OCP()

- Input: None
- Default: None.
- Resets the overcurrent protection (only necessary when the chopping is disabled).


## Overcurrent protection

The CryoSwitch Controller features two overcurrent protection behaviors. The first one occurs when chopping is enabled, chopping essentially clips the current to the maximum set value by “chopping off” the excess current. The second behavior is when the chopping is disabled, in this case, as soon as the overcurrent condition is met, the output gets disabled in a latching way.

The following image shows the current waveforms when the chopping function is enabled and disabled.

The **blue** trace shows the case when the chopping feature is enabled. As soon as the current exceeds the preset threshold (60mA in this case), the controller reduces the voltage in order to reduce the current, also called “Chopping”. After the current drops, the controller will increase the voltage in order to increase the current. This cycle continues until the end of the pulse.

On the other hand, the **orange** trace shows the case when the chopping feature is disabled, as soon as the current exceeds the preset threshold the output is disabled.

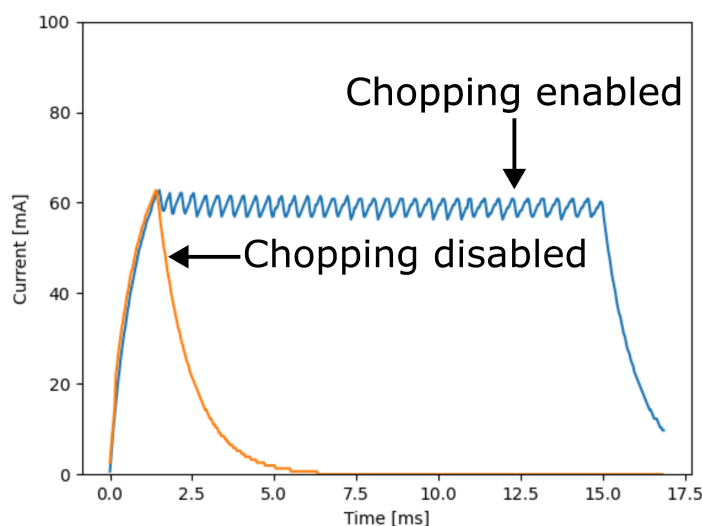

**Figure 7: OCP modes of operation**

## Other functions

### - **set\_sub\_net\_mask(mask='255.255.255.0')**

- Input:
  - Subnet mask number as a string and separated by dots.
- Default: '255.255.255.0'.
- Sets the controller Subnet Mask.

### - **get\_sub\_net\_mask()**

- Input: None
- Returns: Subnet mask number as a string and separated by dots.
- Queries the current controller Subnet Mask.

### - **set\_ip(add='192.168.1.101')**

- Input:
  - IP address number as a string and separated by dots.
- Default: '192.168.1.101'.
- Sets the controller IP address.

### - **get\_ip()**

- Input: None
- Returns: IP address number as a string and separated by dots.
- Queries the current controller IP address.

### - **get\_pulse\_history(port=None, pulse\_number=None)**

- Inputs:
  - *port*: string, optional (A, B, C, D). Filter by the specified port.
  - *pulse\_number*: int, optional. Number of pulses to display.
- Default: None.
- Prints the pulse history, if no pulse\_number is specified it'll show the last 5 pulses.
- Returns: None.

### - **get\_switches\_state(port=None)**

- Inputs:
  - *port*: string (A, B, C, D). Filter by the specified port.
- Default: None.
- Returns the last known state of the port based on the tracking.
- Returns: Dictionary.

**- select\_switch\_model(model='R583423141')**

- Input:
  - Switch the model as a string. Either 'R583423141' or 'R573423600'
- Returns:
  - Bool: True if the selected model is valid
  - Bool: False otherwise.
- Selects the switch model to be used.

**- get\_power\_status()**

- Input: None
- Returns: Power status as an integer, 1 enabled; 0 disabled.
- Queries the converter power status

**- disconnect\_all(port)**

- Inputs:
  - *port*: string (A, B, C, D).
- Default: None.
- Disconnects all the contacts of the specified port.
- Returns: None.

**- smart\_connect(port, contact)**

- Input:
  - *port*: string (A, B, C, D).
  - *contact*: (1,...,6)
- Default: None.
- Disconnects the previously connected contacts and connects the specified switch contact. Based on the tracking history.
- Returns: Current waveform of the specified contact.

**- flash()**

- Input: None
- Returns: None
- Loads the latest firmware to the controller